



**office and laboratory
supplies**

office and laboratory supplies



Our **LEOCRIL**® Arcobaleno range of products is composed of an outbreak of colours to be mixed each other, fluorescent additive and glitters, accessories and **a series of coloured motifs to be inserted in the appliance** which will make the little patients more motivated.




- **diversified particle size of the powder for easy mixing**

- **low residual of free monomer**
to get bubble free polymerized acrylic plates

- **prolonged processing time**
for careful works without any exothermic reaction

- **extremely reduced shrinkage of the two components**
accurate fit of the prostheses in the mouth

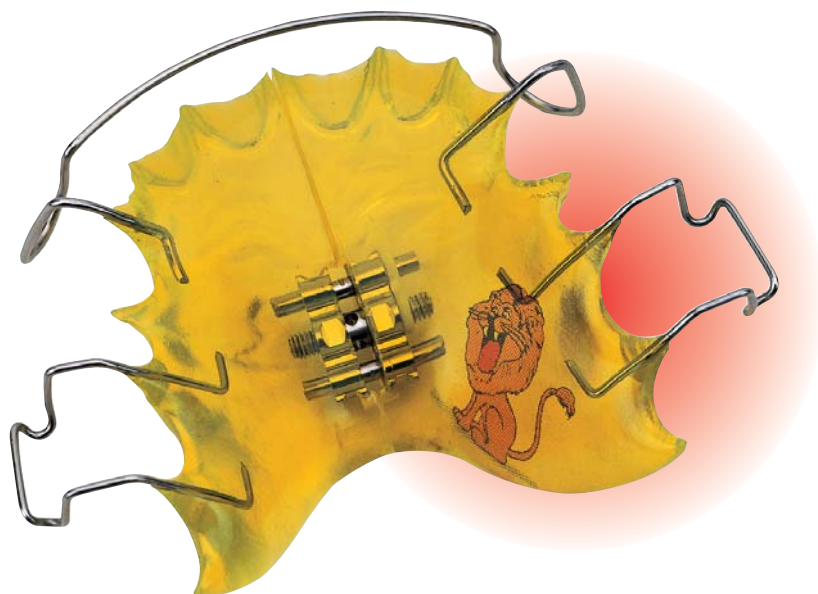

- **utmost precision**
excellent friction of the expansion screws is guaranteed

- **special hardness**
for easy polished surfaces



LEOCRIL[®] orthodontic self-curing acrylic resin

Specific for orthodontic use, either for spray-on or doughing techniques. The special composition of liquid monomer has careful processing properties, without priming of the curing process (exothermic reaction) which may cause an undesired evaporation of monomer. The use of a Leone hydroflask or pressure pot for polymerization at 40° C with a pressure of 2,5 Atm for 20 minutes is suggested. Hard and perfectly finished plates are obtained as no otherwise possible with other methods.

Regular kit

- R6222-95** A - clear
A - clear liquid
clear powder
- R6223-95** B - transparent pink
B - transp. pink liquid
clear powder

The kit includes:
1 Kg powder, 0,5 lt liquid, spray flasks for powder/liquid, funnels, beakers and insulator.

Introductory kit

- R6202-93** A - clear
A - clear liquid
clear powder
- R6203-93** B - transparent pink
B - transp. pink liquid
clear powder

Necessary in the dental office.

It includes: 60 g powder and 45 ml liquid.



LEOCRIL[®] liquid monomer

The package of 6 lt is composed of 3 bottles of 2 lt each for better handling of product

colour		pkg.	
clear	A	0,5 lt	R6122-00
		6 lt	R6142-01
pink transp.	B	0,5 lt	R6123-00
		6 lt	R6143-01



LEOCRIL[®] powder polymer

Clear colour.

R6032-00 Package of 1 Kg

R6042-00 Package of 10 Kg

acrylic resin



LEOCRYL® "Arcobaleno" orthodontic self-curing resin

"Arcobaleno" **LEOCRYL**® method allows the construction of coloured and personalized appliances for a better psychological motivation of the patient.

It includes:

- 1 colour guide with ten coloured plates, 1 fluorescent, 1 glitters and 1 transparent. Each plate is marked with a letter for easy choice of the liquid
- ready-to-use **LEOCRYL**® coloured liquids each marked by a letter
- 1 range of 28 funny motifs for use with the acrylic plate
- 1 fluorescent additive which makes the appliance visible in the dark
- coloured glitters giving a "metallic" effect to the plates.

LEOCRYL® "Arcobaleno" - liquid

			pkg.		pkg.
red	D	R6124-00	250 ml	R6134-00	0,5 lt
purple	E	R6125-00		R6135-00	
magenta	F	R6126-00		R6136-00	
orange	G	R6127-00		R6137-00	
fluorescent yellow	H	R6128-00		R6138-00	
green	I	R6129-00		R6139-00	
turquoise	L	R6130-00		R6140-00	
blue	M	R6131-00		R6141-00	

R6000-93

LEOCRYL®

"Arcobaleno" assorted kit

It includes: 2 Kg clear Leocryl powder, 8 coloured liquids, 8 spray flasks for liquid, 5 packs of glitters, 1 pack of fluorescent additive, funnels, 1 coloured motifs for plates kit and 10 colour guides.



R6000-00

LEOCRYL®

"Arcobaleno" colour guide

Colours shown are only indicative.
Colours may slightly change compared to those obtained with **LEOCRYL**® acrylic resin.

acrylic resin accessories



R4401-93 Coloured motif for plates assorted kit

It includes 28 different motifs, one positioning tweezers and illustrated pamphlet for patient choice.

Assorted kit of 280

Refills of motifs

Packs of 10



Glitters

- R4300-00 Silver
- R4301-00 Gold
- R4302-00 Multicolor "Z"
- R4303-00 Blue
- R4304-00 Red

Nice looking orthodontic plates are obtained by mixing a small quantity of glitters with acrylic powder. Glitters make plates X-ray visible.

Packs of 50 g and spray flask

acrylic resin accessories



R4310-00
"X" fluorescent additive

A fluorescent effect is given by mixing just 15% approx. of product with polymer. The best result is achieved with the use of clear powder and liquid.

Pack of 100 g



R6320-00
Insulator liquid for resin

Pack of 1 lt



R4550-00
Polishing fluid for resin

Spread some drops on the polishing brush.

Pack of 0,5 lt



Accessories

R4401-00 Spray flask for powder

R4405-00 Funnel for liquid

R4409-00 Spray needle for liquid

R4402-00 Spray flask for liquid

R4406-00 Funnel for powder

R4410-00 Arrows

R4403-00 Large beaker

R4407-00 Brush

R4404-00 Small beaker

R4408-00 Spray nozzle for powder/liquid

plaster for orthodontic models and accessories



Extra hard plaster

R5005-00 Bucket of 5 Kg

R5025-00 Bucket of 25 Kg


Orthodontic plaster type III, snowy white colour.

Mixing ratio: 26-28 ml: 100 g.

Fine detailed, smooth and hard models (over 650 Kg/cm²) are obtained.

Easy to polish with R4520-00.

Ideal for making study models and plaster casts for collection.

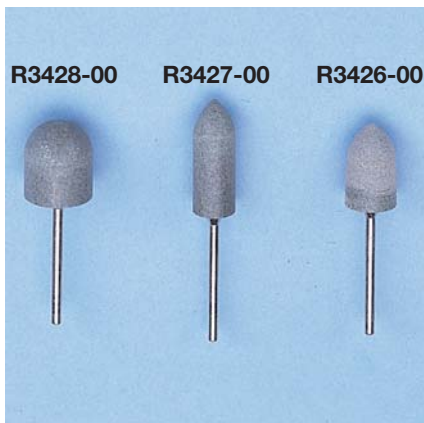

Silicon mounted polishers

For perfect polishing of plates. They do not damage hooks or arch wires. All with 2,35 mm diameter shaft.

Packs of 6




Gloss for plaster models

R4520-00 Pack of 0,5 lt

R4530-00 Pack of 5 lt

It is necessary to use a good ortho plaster but also a good final model gloss to obtain a well-made study model. This model gloss provides polish surfaces with fine details and protects them over the time.

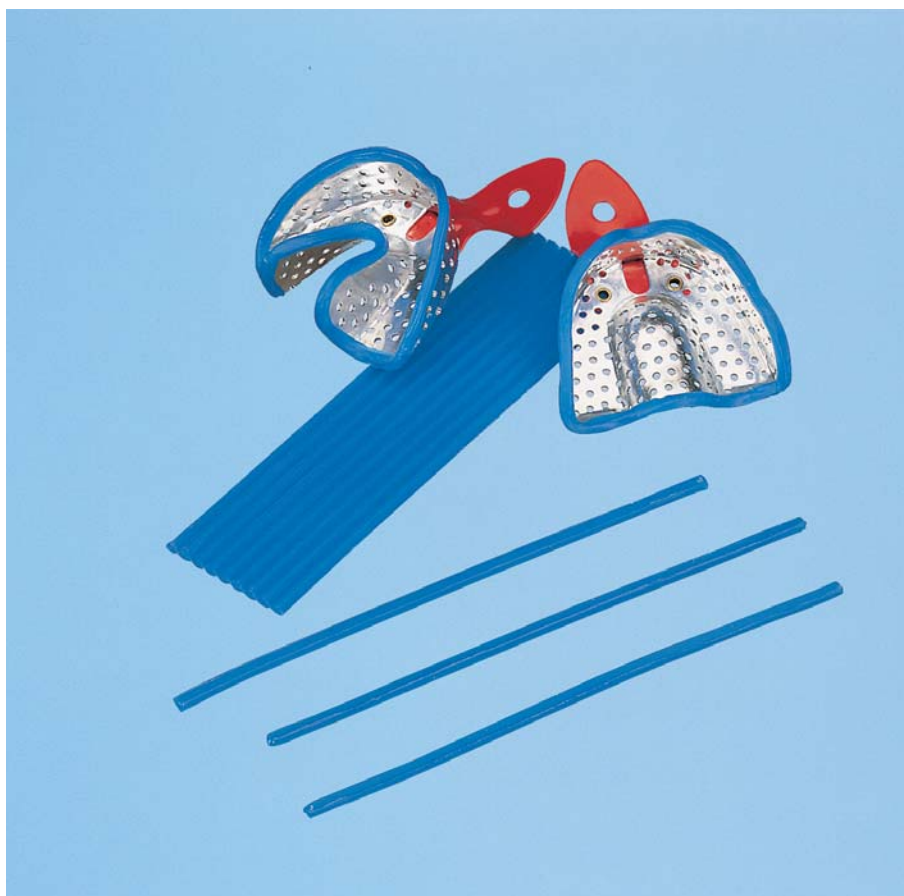

R3712-00

Quadrangular wax sticks

Wax in 30 cm long sticks, 3x3 mm cross section.

The special adherence and high workability make it useful to border the impression trays to respect the edges of the functional impression.

Pack of 90 sticks



R3710-00

Utility wax

Highly soft wax in strips which bends at room temperature without breaking. Thanks to its good adherence and high workability, it serves to border the impressions to obtain the base of the mould.

Pack of 24 strips



R3711-00

Modelling wax

Wax in sheets. It moulds, both cold and hot, without breaking.

Pack of 500 g



R3713-00

Sticky wax

High sticky power. Thanks to its extreme strength of adhesion, it allows fixing of the various elements (screws, wires etc.) to the model.

Pack of 100 g in stick



R3730-00

Bite wax

Necessary for taking accurate occlusal impressions.
Specially aluminium foil laminated, it lessens distortions.
Pleasantly scented, 6 mm thick.

Pack of 24



R3714-00

Thermowax for undercuts

High melting point: 104° C. Particularly suitable to eliminate undercuts.
It prevents damage to be caused to the acrylic brightness.

Pack of 100 g in stick



R3716-00

Natural protective wax

Specially intended for orthodontic use.
The natural composition of this wax makes it to be absolutely safe and if swallowed, it does not cause any illness. Mint flavoured.

Pack of 10 boxes, 5 sticks each

silver solders and accessories



Solders for stainless steel

R0220-00 Plate 0,43 mm thick

R0221-00 Wire form, 1 mm in diameter

For use with bridges and crowns.
Cadmium free. Silver content: 27%.
Melting point: 750° C.

Packs of 10 g approx.



Orthodontic solders

R0222-00 Wire form, 0,40 mm in diameter

R0223-00 Tape form, 1 x 0,5 mm

R0224-00 Rods, 1 mm in diameter

For all orthodontic brazing works.
Cadmium free. Silver content: 55%.
Low melting point: 600° C.

Packs of 10 g approx.



R0225-00 Solder powder

Powder for solders R0220-00 and R0221-00.

Pack of 50 g approx.



R0226-00 Flux

In paste, for solders:
R0222-00, R0223-00 and R0224-00.

Pack of 30 g approx.



R0227-01 Thermo insulator paste for solders

It protects plaster and resin during soldering procedure. Excellent consistency. Easy to remove with water.

Pack of 235 ml